

“FINALLY A
WHEELCHAIR
THAT GETS ME
TO THE TOP OF
MY DRIVEWAY.
SOMETIMES
IT’S THE LITTLE
THINGS LIKE
BEING ABLE TO
COLLECT MY OWN
MAIL THAT MAKE
A DIFFERENCE”

Sue
United Kingdom

FULL LIGHTING PACKAGES

Been out having too much
fun and it’s getting late?
Put your lights on.

TAKE CONTROL

An extensive range of
advanced control systems
and environmental controls
are available. Easy to
program and use.

REHAB SEATING

Infinite options for infinite
body shapes. Our rehab
seating range is easily
configurable to suit your
individual needs.

TAILORED FOR YOU

Select the right headrest,
backrest, armrests, legrests,
seating and cushioning to
suit your lifestyle.

magic
Mobility®

Frontier V4 FWD
**FOR THE
ROAD LESS
TRAVELLED**



CE UK
CA

V4/AUST/AUGUST26



Frontier
V4®



FRONTIER V4 FRONT WHEEL DRIVE

A powered wheelchair for people who know what they want. The V4 FWD offers a great solution to some of life's most tricky corners. The thinking person's chair, the V4 FWD is designed for those wanting to get around their homes with ease, using desks and tables, or driving frequently. Choose Off-road, Crossover or Urban wheels to ensure you get the best powerchair for your lifestyle.

Basic specifications:*

- Length: 995mm
- Width: 710mm - AT wheels
700mm - Crossover wheels
640mm - Urban wheels
- Seat to floor height: 435mm
- Base turning radius: 520mm
- Weight limit: 182kg (400lbs),
160kg (350lbs) with lift and/or tilt
- Top speed: 10kph (6.2mph)
- Motors: 800W Gear in Line
- Wheels: Off-road 14" black,
Crossover 14" black or grey,
Urban 14"x3" black or grey
(solid or pneumatic)
- Castors: 220mm (2.80/2.50-4),
black or grey (solid or pneumatic)
- Gel batteries
- Penny and Giles controller with
LCD display (R-Net)
- Centre post legrest with one piece
footplate
- Rehab seat pan (from 12"-24" wide)
- Crash test compliant tie down
points

Rehab positioning:

- Power options: 300mm (12")
seat lift, tilt (45 degrees), lift and
tilt, transfer lift and tilt only, transfer
tilt with lift and tilt, power recline
and more
- Seating and backrests: MPS, rehab
backrest, adjustable seat pans and
seating
- Headrests and mounts: MPS,
rehab
- Legrests: centre post, side mount,
swingaway, power, gas or manual
elevating
- Footplates: 1 or 2 piece fixed/flip
up/angle adjustable, any specific
requirement
- An extensive range of armpads,
supports, restraints and pads are
available

The right wheelchair for:**

- Adventurous off-road recreation
- Country living and working
- Urban outdoor living and working
- Roomy indoor spaces

Options and accessories:

- One-click activator
(CJSM 2 or JSM joystick)
- Fold forward back
- Retractable docking pin
- Puncture resistant tyre liners
(Off-road and Crossover
wheels only)
- 24V LED lighting kit
- 12V accessory socket and
2 x USB charger
- Wheelchair jack
- Tray tables
- Luggage rack, bag hooks, bags,
pockets
- Special controls: attendant,
head, chin, foot, sip and puff,
environmental controls
- RAM® mounting system (tablet
and swing away tablet holder,
mobile phone holder, swing away
camera mount, attendant cup
holder, action camera holder and
cup holder)

Base and tyre options:



SPACE BLUE BASE
BLACK OFF-ROAD TYRES & RIMS



LUCENT LIME BASE
BLACK CROSSOVER TYRES



BURNT ORANGE BASE
BLACK URBAN TYRES

"WE'VE BEEN
CREATING
HIGH-QUALITY
WHEELCHAIRS
FOR 30 YEARS,
SO THAT YOU
CAN LIVE THE
LIFE YOU
CHOOSE."

Ashley Daff
R&D Director




**ATTENTION
TO DETAIL**

We know that small touches
make a big difference.
Our powerchair range is
configurable and modifiable
to suit your individual needs.

**ROBUST
QUALITY**

We never compromise on
quality. Your powerchair
is your mobility, therefore
reliability is paramount.

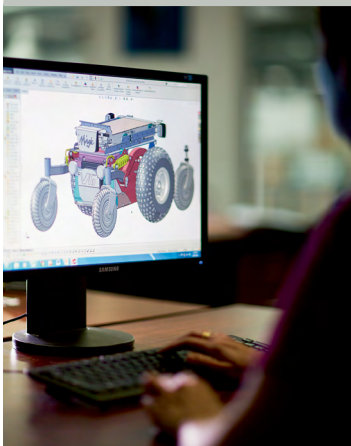

**ENGINEERED TO
THE LAST BOLT**

We listen to and value
your feedback. Our team
of in-house engineers work
side by side with end users
and therapists. We're all
about innovation; if we
think we can make our
powerchairs better, we will.

**AUSTRALIAN
DESIGNED
AND MADE**

We are proud that our
powerchair range is
designed and made in
our Australian
manufacturing facility.

*Specifications subject to change. Check online for current version
**Please remember that your wheelchair has electric motors and must
not be driven through water, into rivers, creeks or the sea



Magic Mobility +61 3 8791 5600
enquiries@magicmobility.com.au
Explore with us magicmobility.com.au



Scan for agent locations